


Unison UK-C18
 Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18

250 x 3 mm

Application

Sialylglycopeptide in egg yolk extract

卵黄抽出物中のシアリルグリコペプチド

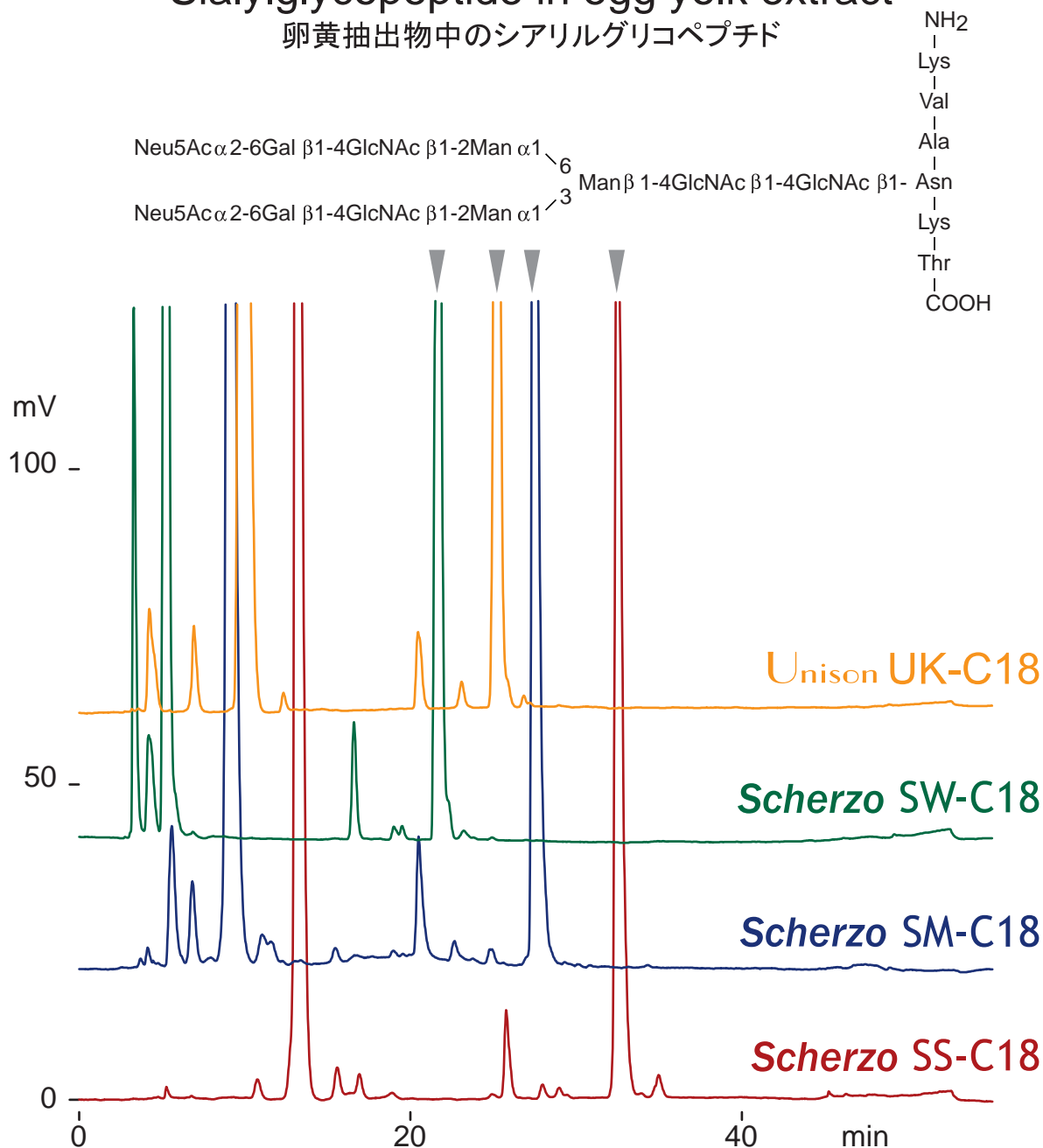

250 x 3 mm

A: 20mM ammonium acetate

B: 100mM ammonium acetate /acetonitrile = 50 /50

0-10 %B(0-40min), 10-100 %B(40-50min)

0.4 mL/min (12MPa), 37 deg.C

ELSD (spray chamber 40 deg.C, drift tube 100 deg.C)

5 μ L

Courtesy of FUSHIMI Pharmaceutical Co., LTD.

Distributor

Imtakt USA (info@imtaktusa.com)

Imtakt Corp./JAPAN (info@imtakt.com)

LCC Engineering & Trading GmbH, Steinbruchstr. 4, CH-4622 Egerkingen, Switzerland.

www.chromatographyshop.com, info@chromatographyshop.com,

Tel: +41 62 398 52 71, Fax: +41 62 398 52 74