

Conditions



Column: ACE Excel 2 C18-PFP
Dimensions: 100 x 2.1 mm
Part Number: EXL-1010-1002U
Mobile Phase: A: 20 mM HFBA in H₂O/MeCN (98:2 v/v)
 B: 20 mM HFBA in MeCN/H₂O (98:2 v/v)
Gradient:

Time (mins)	%B	Curve
0.0	5.0	-
2.0	15.0	6
4.5	19.0	6
5.5	19.5	8
6.0	22.0	6
7.0	35.0	6
9.0	48.0	8
9.5	5.0	6

Flow Rate: 0.4 mL/min
Temperature: 40 °C
Detection: Positive ESI MRM (transitions as shown)
Sample: Extraction at low pH, clean up with WCX SPE cartridge
 Egg sample spiked at 100 µg/kg (CCα)

Analytes

GNT Gentamicin
 TBM Tobramycin
 KNM Kanamycin
 DHS Dihydrostreptomycin
 STR Streptomycin
 HGB Higromycin-B
 SPT Spectinomycin

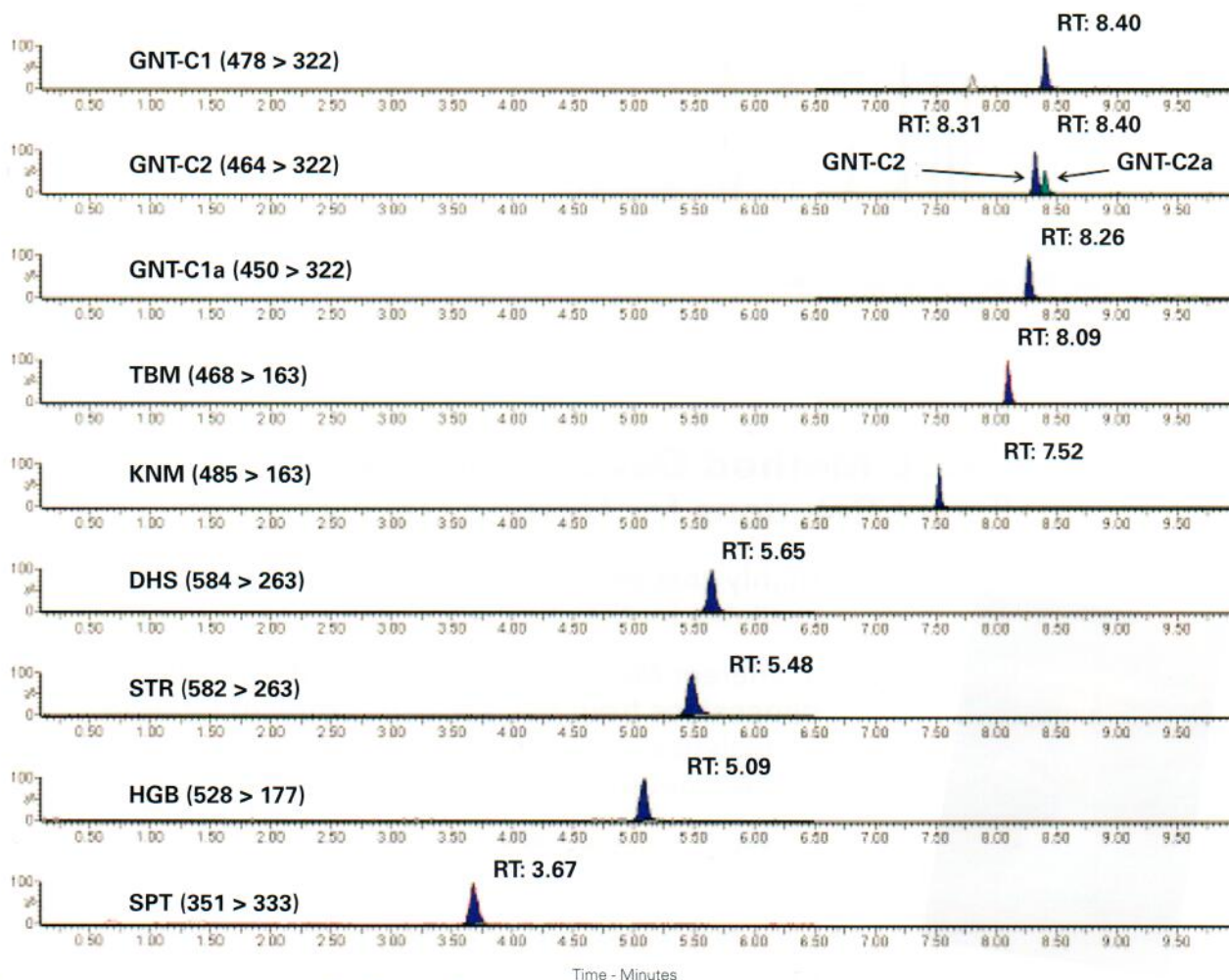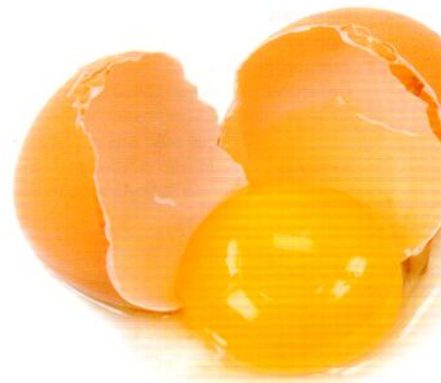

Reproduced with permission of Veterinary Drug Residues Department, Centro Nacional de Alimentación (CNA, AECOSAN), Madrid, Spain

Available from



LCC Engineering & Trading GmbH
 Steinbruchstrasse 4, CH-4622 Egerkingen, Switzerland
 Tel.: +41 62 398 5271, Fax.: +41 62 398 5274
 info@chemsupply.ch
 www.chromatographyshop.com