





Scherzo SM-C18

10 x 3 mm

Application

LC-MS: Adenosine phosphates (AMP, ADP, ATP)

アデノシンリン酸 (AMP, ADP, ATP)

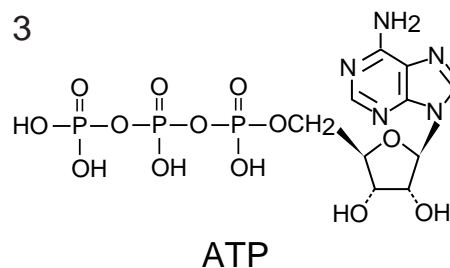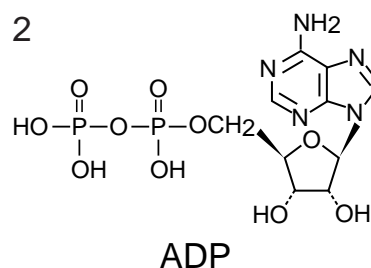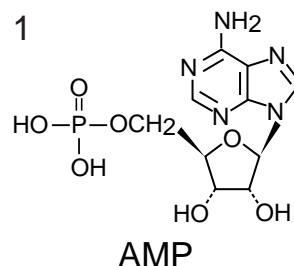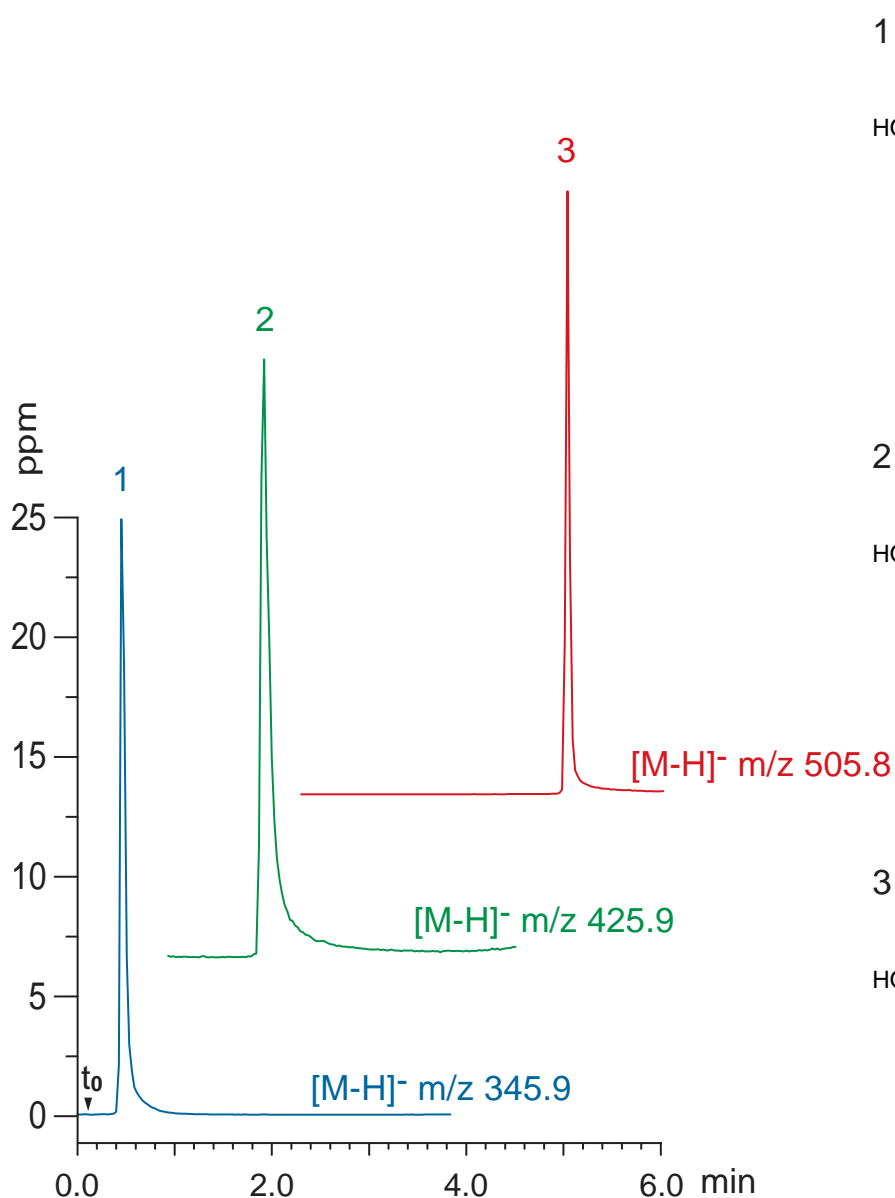

Scherzo SM-C18, 10 x 3 mm

A: water / formic acid = 100 / 0.5

B: (100mM HCOONH₄ : 100mM NH₄OH = 97:3) / MeOH = 90 / 10

0-100 %B (0-2.5 min), 100 %B (2.5-6 min)

0.4 mL/min (1.6 MPa), 37°C, 20 uL (25 ppm, water)

Single Quad. MS (ESI, negative)

Distributor

LCC Engineering & Trading GmbH, Steinbruchstr. 4, CH-4622 Egerkingen, Switzerland.

www.chromatographyshop.com , info@chromatographyshop.com,

n s

Tel: +41 62 398 52 71 , Fax: +41 62 398 52 74