



Scherzo SS-C18  
 Scherzo SM-C18  
 Scherzo SW-C18  
 Unison UK-C18

150 x 3 mm

Application

**Arsenobetaine**

アルセノベタイン

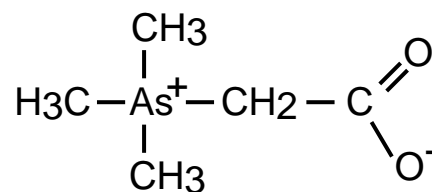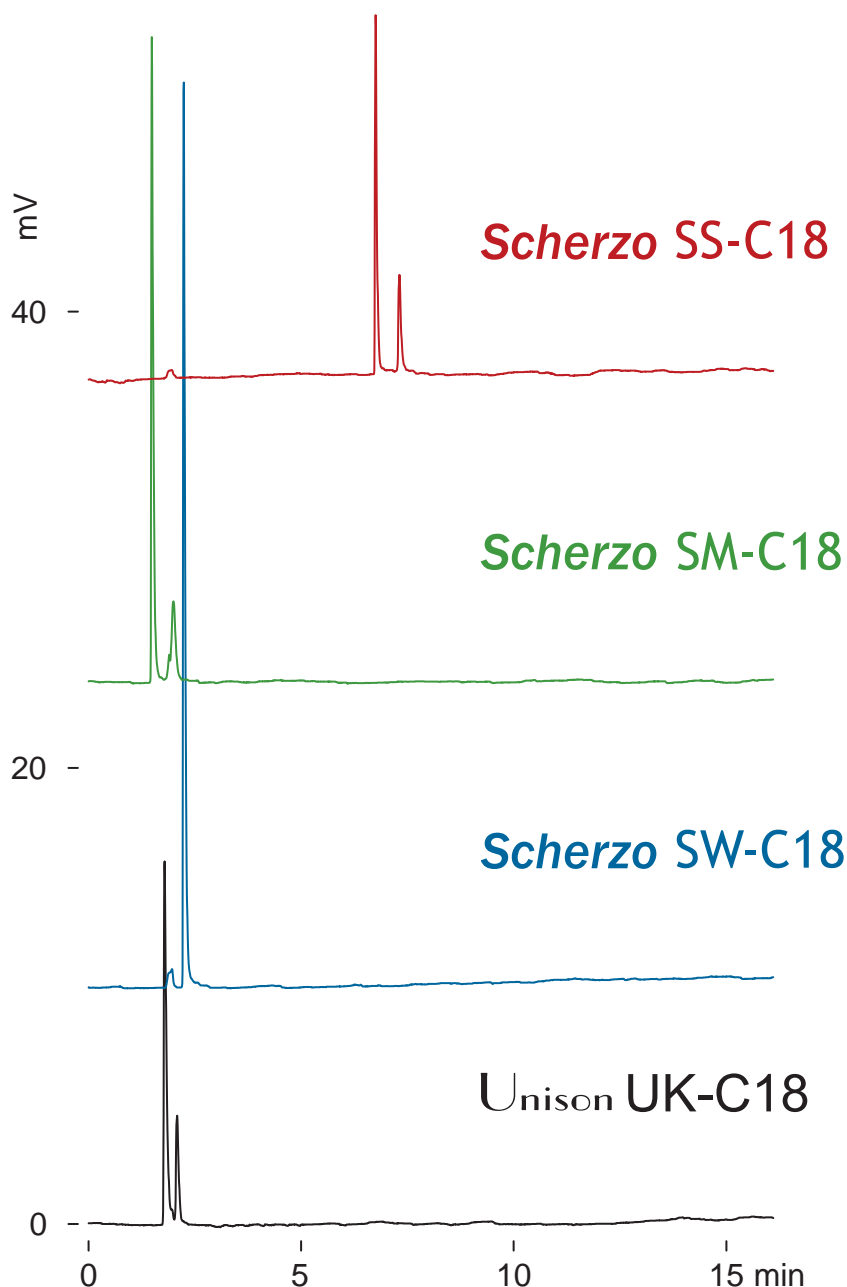

arsenobetaine

150 x 3 mm

A: water / formic acid = 100 / 0.2

B: 100 mM ammonium formate / acetonitrile = 50 / 50

5-80 %B (0-15 min )

0.4 mL/min (8-9 MPa), 37 deg.C, ELSD (spray chamber 50 deg.C, drift tube 100 deg.C)

20 uL ( 0.49 ug )

**Distributor**

LCC Engineering & Trading GmbH, Steinbruchstr. 4, CH-4622 Egerkingen, Switzerland.

www.chromatographyshop.com , info@chromatographyshop.com,

Tel: +41 62 398 52 71 , Fax: +41 62 398 52 74