





Scherzo SW-C18
Unison UK-C18

250 x 3 mm

Application

Coffee

コーヒー

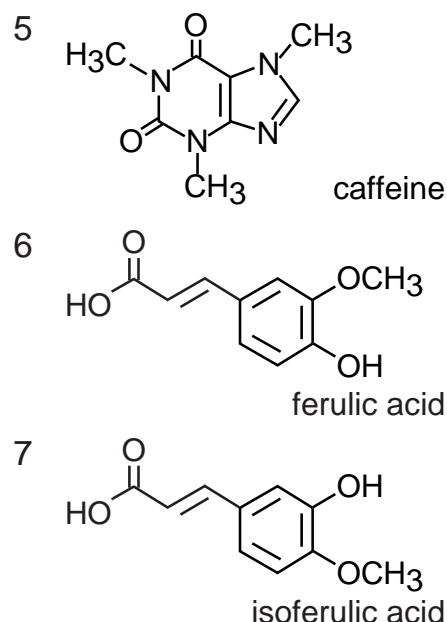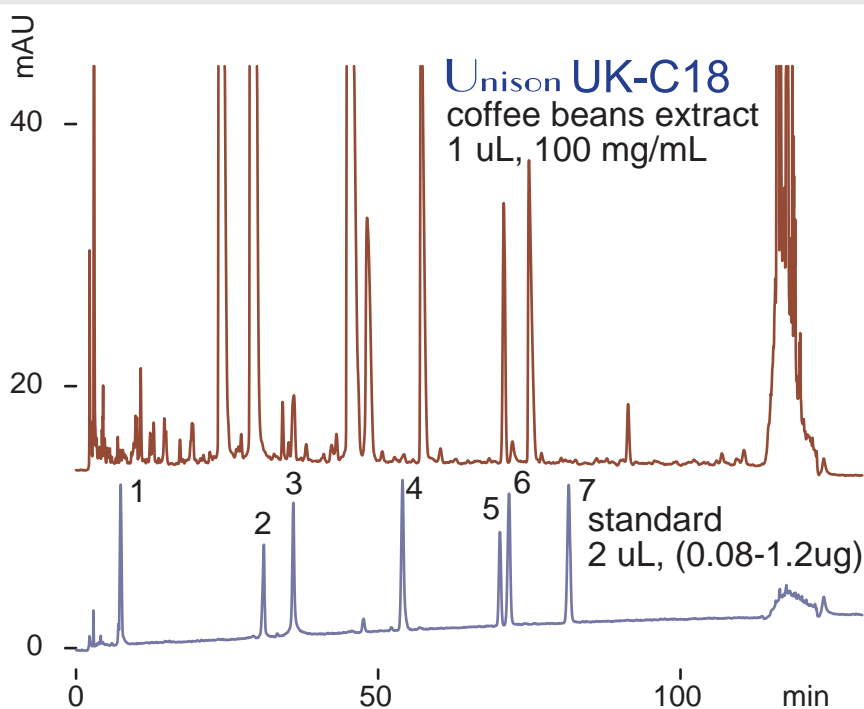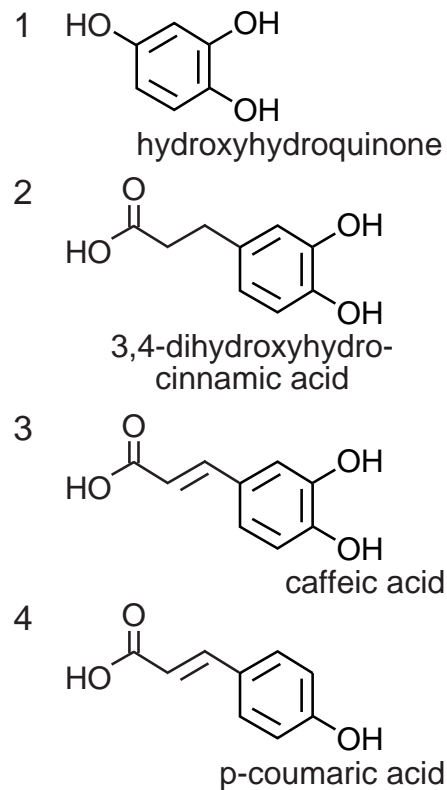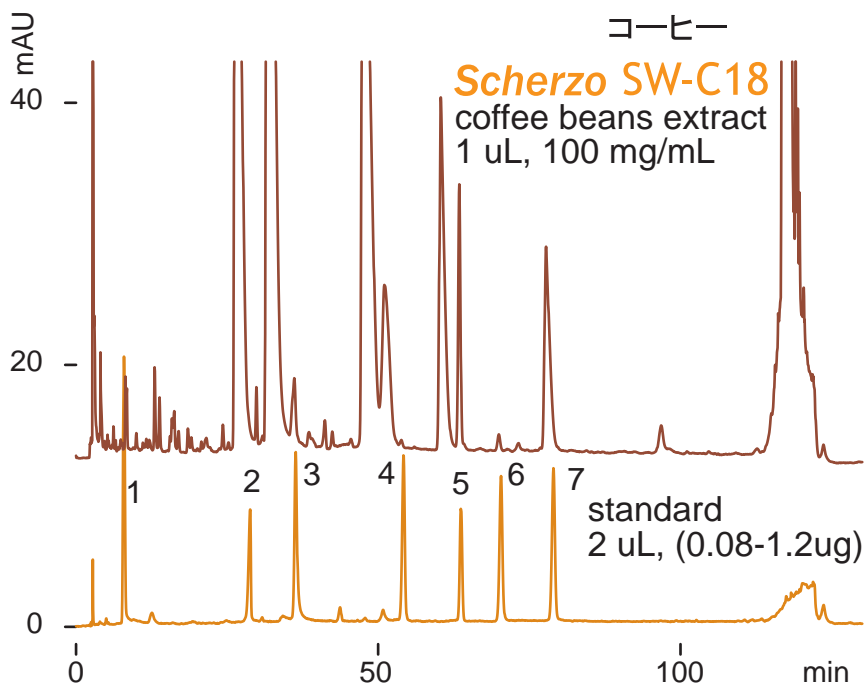

250 x 3 mm

A: 10mM formic acid /10mM ammonium formate = 5 /95

B: (40mM formic acid /40mM ammonium formate = 5 /95) /acetonitrile = 50 /50

0-20 %B (0-110min), 20-100 %B (110-120min)

0.5 mL/min (16MPa), 37 deg.C, 260 nm

Distributor

LCC Engineering & Trading GmbH, Steinbruchstr. 4, CH-4622 Egerkingen, Switzerland.

www.chromatographyshop.com , info@chromatographyshop.com,

Tel: +41 62 398 52 71 , Fax: +41 62 398 52 74