



Cadenza CD-C18
 Unison UK-C18
 Cadenza CL-C18
 Scherzo SM-C18

150 x 3 mm


Application

Endosulfan isomers

エンドスルファン異性体

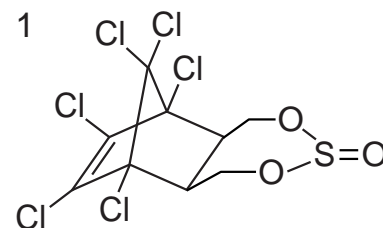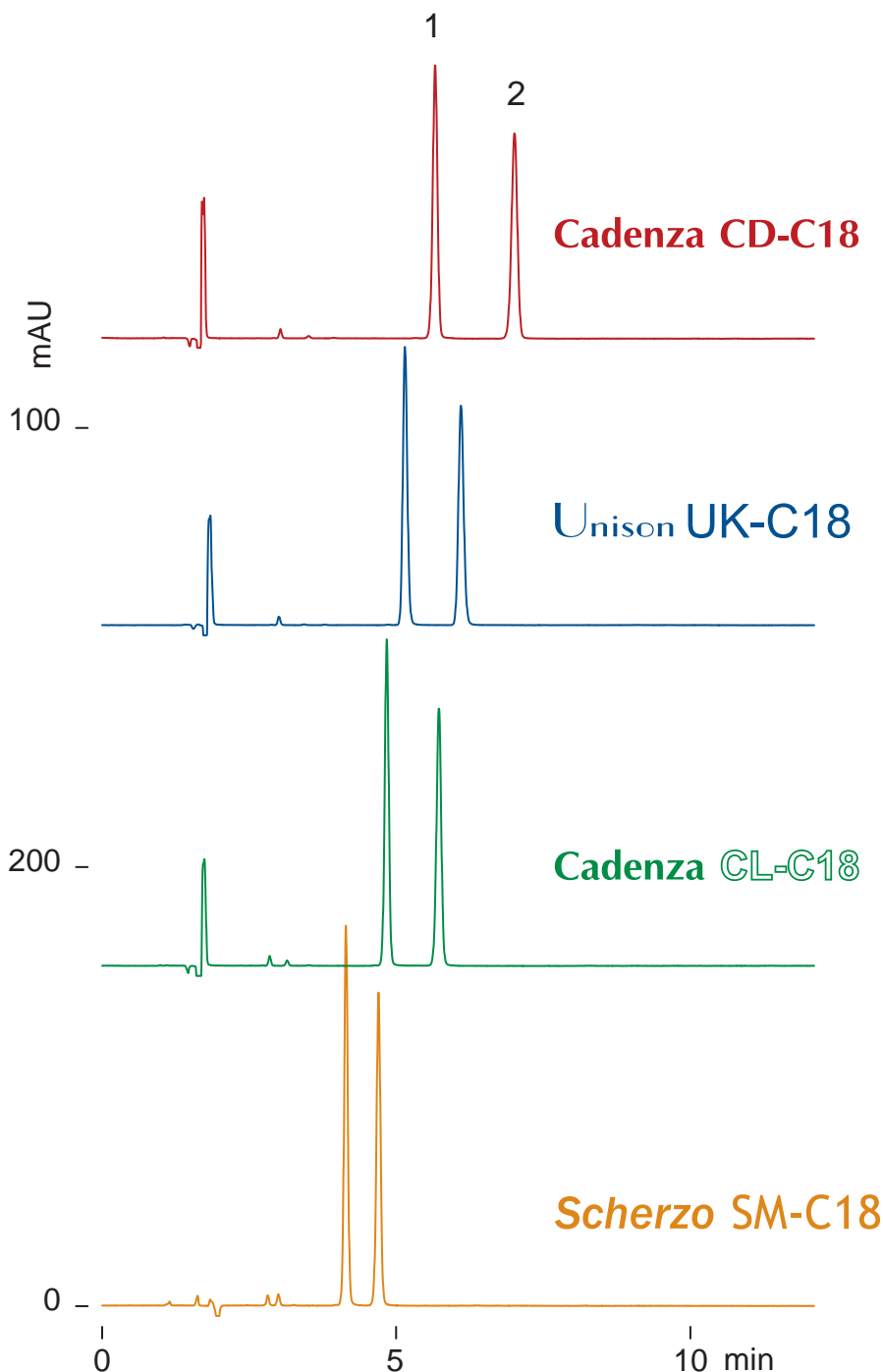

beta-endosulfan

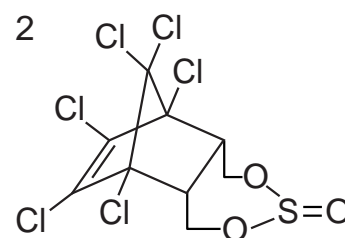

alpha-endosulfan

150 x 3 mm

water / acetonitrile / acetic acid = 20 / 80 / 0.1

0.4 mL/min (5-6 MPa), 37 deg.C, 220 nm, 1.5 uL (0.75 ug, methanol)

Distributor

LCC Engineering & Trading GmbH, Steinbruchstr. 4, CH-4622 Egerkingen, Switzerland.

www.chromatographyshop.com , info@chromatographyshop.com,

Tel: +41 62 398 52 71 , Fax: +41 62 398 52 74