








Scherzo SW-C18

150 x 3 mm

Application

LC-MS: Nicotine related compounds (1)

ニコチン関連化合物 (1)

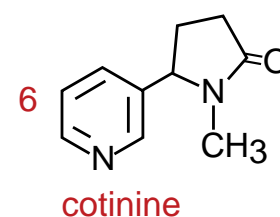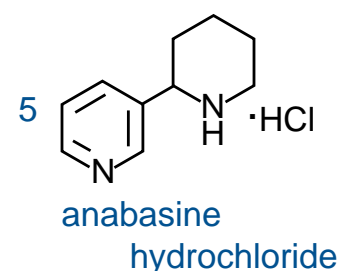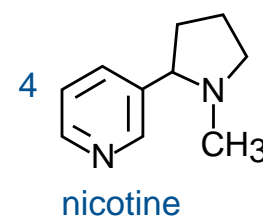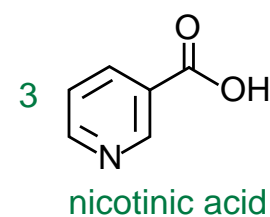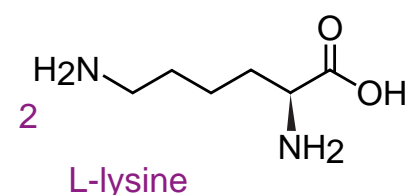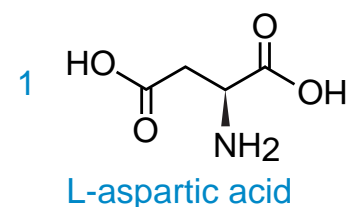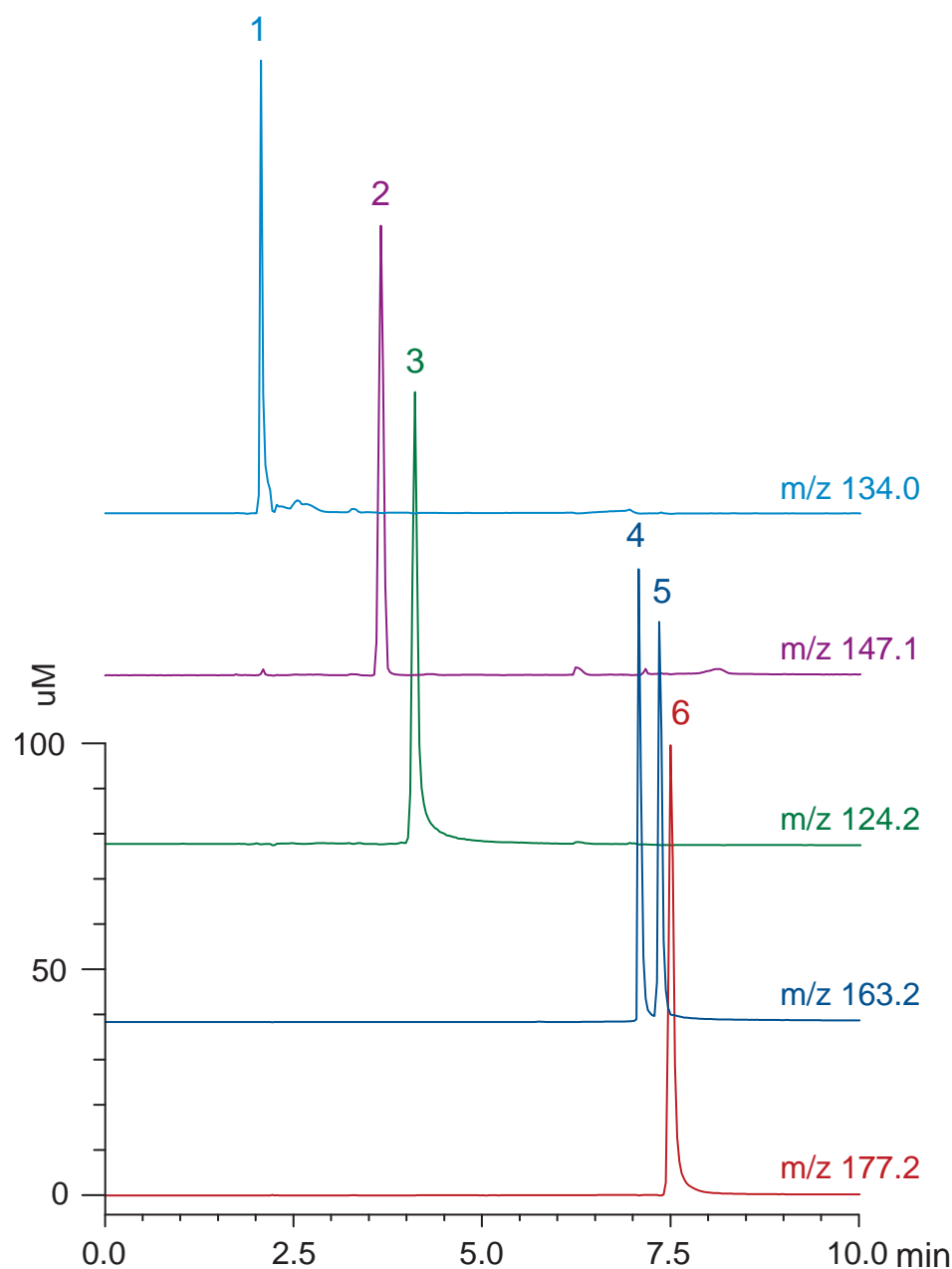

Scherzo SW-C18, 150 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 20 / 80

0-20 %B (0-4 min), 80 %B (4-10 min)

0.4 mL/min (9 MPa), 37°C, 5 uL (100 nmol/mL, 1% HCOOH)

Single Quad. MS (ESI, positive)

Distributor

LCC Engineering & Trading GmbH, Steinbruchstr. 4, CH-4622 Egerkingen, Switzerland.

www.chromatographyshop.com , info@chromatographyshop.com,

Tel: +41 62 398 52 71 , Fax: +41 62 398 52 74