


Scherzo SM-C18

150 x 3 mm

Application

### Nucleotide triphosphates

ヌクレオチド3リン酸

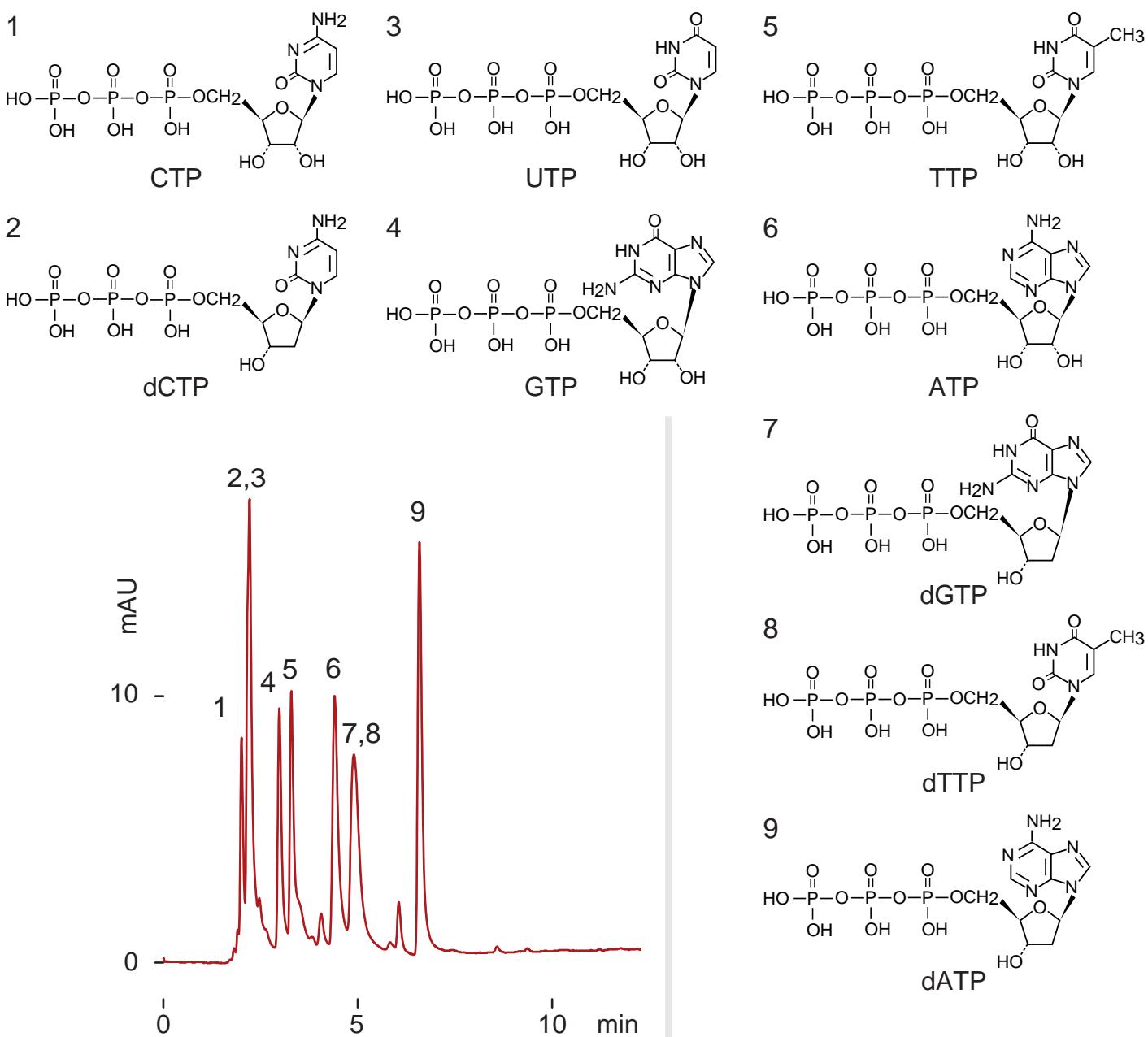

Scherzo SM-C18, 150 x 3 mm

A: 50mM HCOONH<sub>4</sub> : 50mM NH<sub>4</sub>OH = 9 : 1, pH 8.6

B: (100mM HCOONH<sub>4</sub> : 100mM NH<sub>4</sub>OH = 9 : 1, pH 8.6) / acetonitrile = 80 / 20  
0-50 %B (0-10min)

0.4 mL/min (8MPa), 37 deg.C, 260 nm, 1 uL (0.1ug, water)

#### Distributor

LCC Engineering & Trading GmbH, Steinbruchstr. 4, CH-4622 Egerkingen, Switzerland.

[www.chromatographyshop.com](http://www.chromatographyshop.com), [info@chromatographyshop.com](mailto:info@chromatographyshop.com),

Tel: +41 62 398 52 71, Fax: +41 62 398 52 74