

Sugars – Maple Syrup


Application #AN4310

Conditions

Column: ACE Excel 1.7 NH₂
Dimensions: 100 x 3.0 mm
Part Number: EXL-1714-1003U
Mobile Phase: MeCN/H₂O (70:30 v/v)
Flow Rate: 0.43 mL/min
Injection: 10 µL
Temperature: 35 °C
Detection: RI, 35 °C
Sample: 100 µL maple syrup in 9900 µL mobile phase
System: Chromaster 600

Analytes

1. Fructose
2. Glucose
3. Sucrose

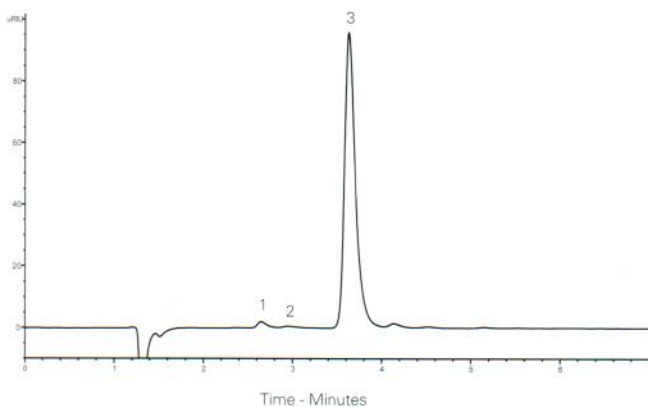

Sugars – Monosaccharides


Application #AN4330

Conditions

Column: ACE Excel 5 NH₂
Dimensions: 150 x 4.6 mm
Part Number: EXL-1214-1546U
Mobile Phase: MeCN/H₂O (75:25 v/v)
Flow Rate: 1 mL/min
Injection: 10 µL
Temperature: 35 °C
Detection: RI, 35 °C
Sample: 5 mg/mL per monosaccharide in MeCN/H₂O (70:30 v/v)
System: Chromaster 600

Analytes

1. Fructose
2. Galactose
3. Rhamnose

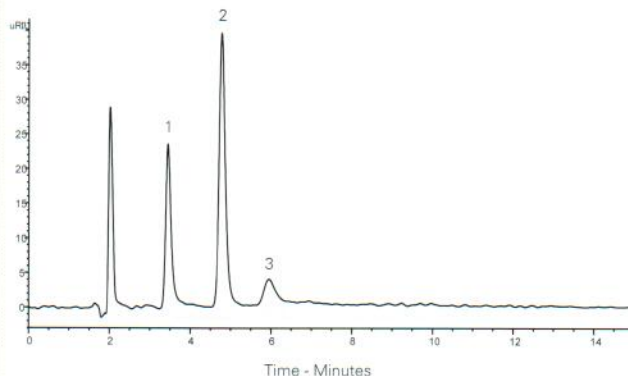

Please contact us for further information and advice on specific applications or for method development support

Available from



LCC Engineering & Trading GmbH

Steinbruchstrasse 4, CH-4622 Egerkingen, Switzerland

Tel.: +41 62 398 5271, Fax.: +41 62 398 5274

info@chemsupply.ch

www.chromatographyshop.com