


Scherzo SM-C18

150 x 2 mm

Application










Analysis of Water-soluble Vitamins using LCMS-2020

LCMS-2020を用いた水溶性ビタミンの分析

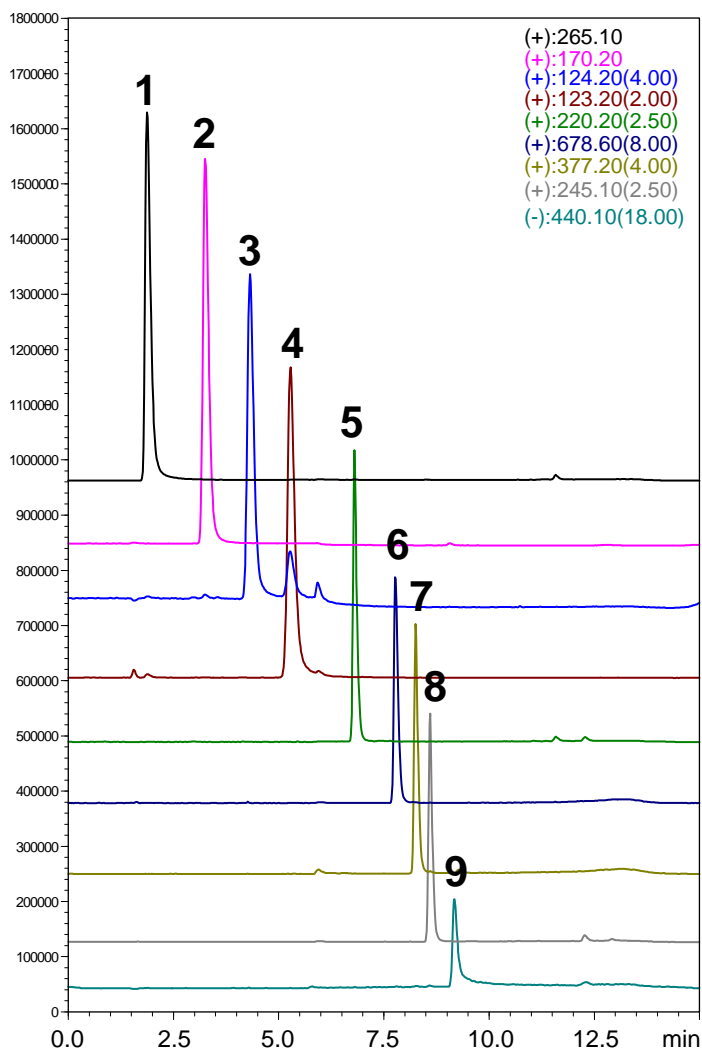

Shimadzu LCMS-2020
 (ESI, positive/negative, SIM)

Scherzo SM-C18, 150 x 2 mm
 A: 5 mmol/L ammonium formate
 + 0.1% formic acid - water
 B: acetonitrile
 0 - 55 %B (0-10 min), 0 %B(10.01-20 min)
 1,000 pg/uL, 1uL

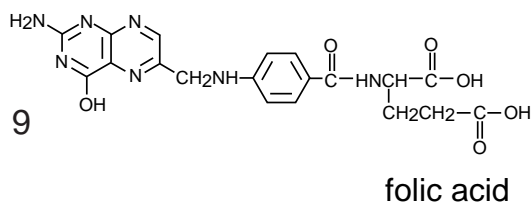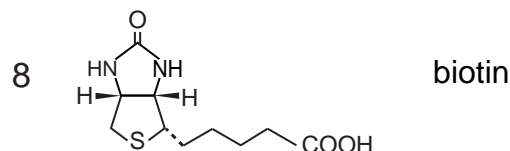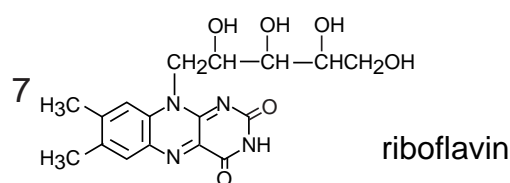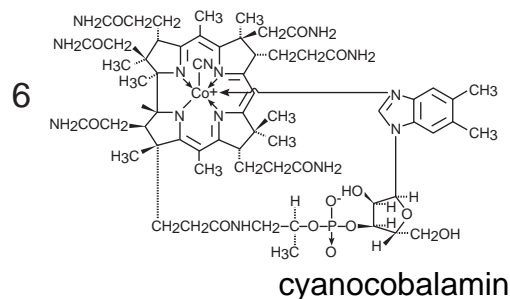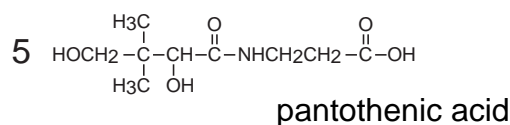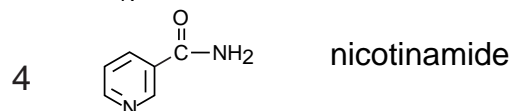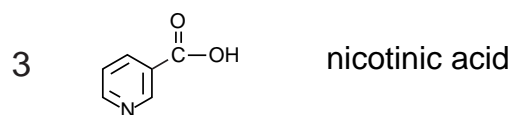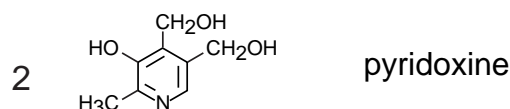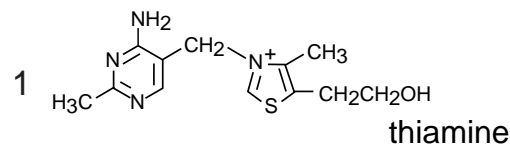

Distributor

LCC Engineering & Trading GmbH, Steinbruchstr. 4, CH-4622 Egerkingen, Switzerland
info@chromatographyshop.com, www.chromatographyshop.com Phone: +41 62 398 52 71, Fax: +41 62 398 52 74