


Conditions


Column: ACE Excel 3 C18-Amide
Dimensions: 250 x 2.1 mm
Part Number: EXL-1112-2502U
Mobile Phase: 40 mM ammonium phosphate pH 2.5 in H₂O
Flow Rate: 0.21 mL/min
Injection: 5 µL
Temperature: 25 °C
Detection: UV, 214 nm

Analytes


1. Oxalic acid
2. Tartaric acid
3. Malic acid
4. Lactic acid
5. Ascorbic acid
6. Citric acid

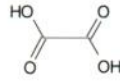


Oxalic acid

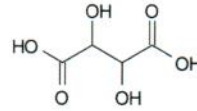


Tartaric acid

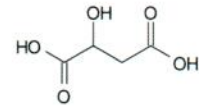


Malic acid

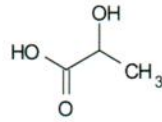


Lactic acid

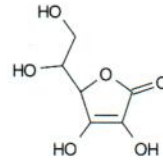

Ascorbic acid

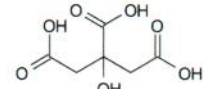

Citric acid

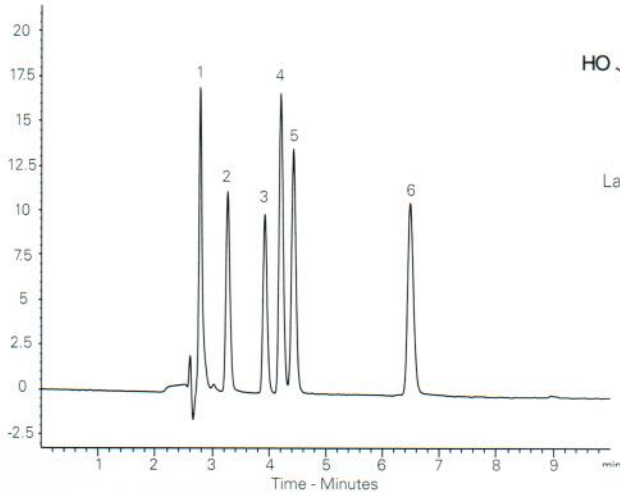

Available from



LCC Engineering & Trading GmbH
 Steinbruchstrasse 4, CH-4622 Egerkingen, Switzerland
 Tel.: +41 62 398 5271, Fax.: +41 62 398 5274
info@chemsupply.ch
www.chromatographyshop.com

