


**Intrada WP-RP**  
Presto FF-C18


250 x 4.6 mm

Application

**Collagen**  
コラーゲン

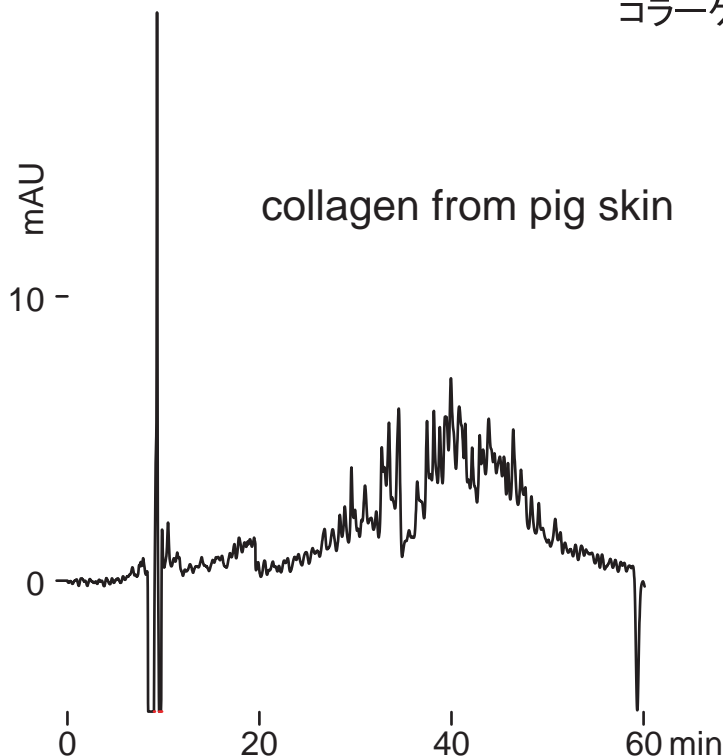

**Intrada WP-RP**, 250 x 4.6 mm  
A: water /trifluoroacetic acid = 100 /0.1  
B: water /acetonitrile /trifluoroacetic acid = 40 /60 /0.1  
15-25 %B (0-50min)  
0.4 mL/min (3MPa), 80 deg.C, 220 nm  
100 uL (0.01%, 1mM HCl)

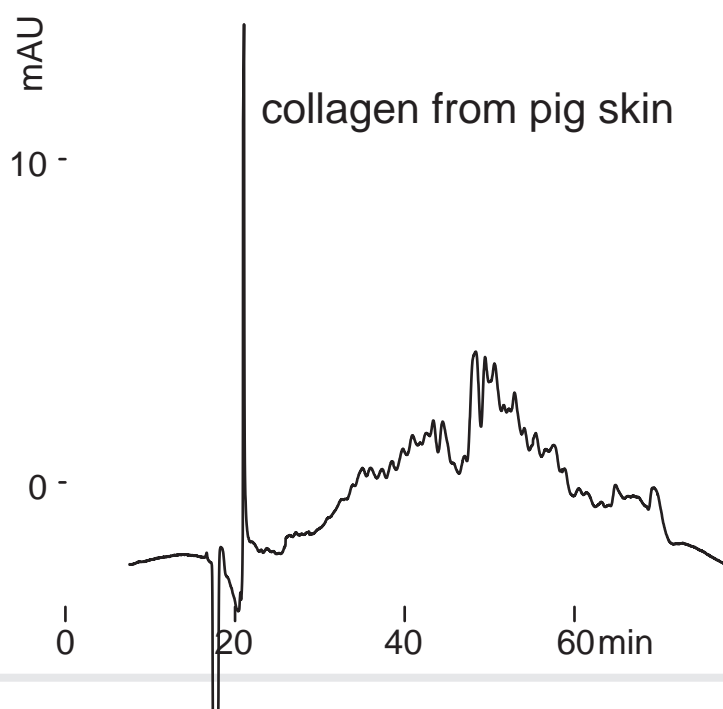

**Presto FF-C18**, 250 x 4.6 mm  
A: water /trifluoroacetic acid = 100 /0.1  
B: water /acetonitrile /trifluoroacetic acid = 40 /60 /0.1  
35-45 %B (0-50min)  
0.15 mL/min (7MPa), 80 deg.C, 220 nm  
50 uL (0.01%, 1mM HCl)

Available from

**LCC Engineering & Trading GmbH**  
Steinbruchstrasse 4, CH-4622 Egerkingen, Switzerland  
Tel.: +41 62 398 5271, Fax.: +41 62 398 5274  
[info@chemsupply.ch](mailto:info@chemsupply.ch)  
[www.chromatographyshop.com](http://www.chromatographyshop.com)

**Intakt Corp./JAPAN** ([info@imtakt.com](mailto:info@imtakt.com))  
Other Countries

TI686E-KH29