


Presto FF-C18

250 x 4.6 mm


Application

Hyaluronic acid

ヒアルロン酸

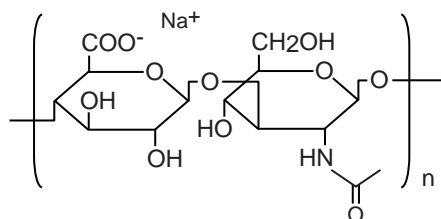

hyaluronic acid sodium salt
(ca. 1 MDa)

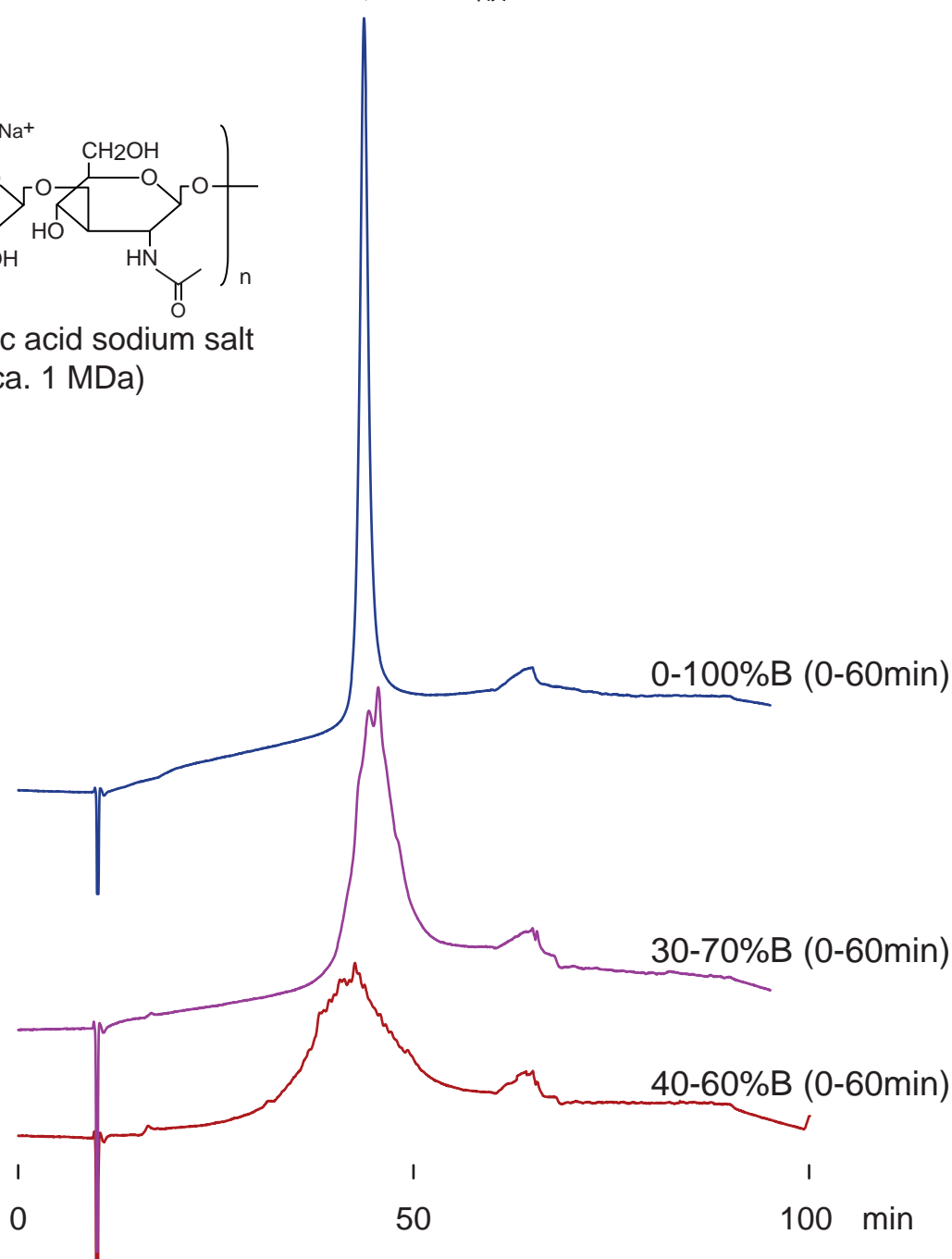

Presto FF-C18, 250 x 4.6 mm

A: water /acetonitrile /trifluoroacetic acid = 95 /5 /0.1

B: water /acetonitrile /trifluoroacetic acid = 90 /10 /0.1

0.15 mL/min (0-60min), 0.3mL/min (60-90min, initialization)

60 deg.C, 210 nm, 4 uL (8ug, 0.05%TFA)

Available from

LCC Engineering & Trading GmbH
Steinbruchstrasse 4, CH-4622 Egerkingen, Switzerland
Tel.: +41 62 398 5271, Fax.: +41 62 398 5274
info@chemsupply.ch
www.chromatographyshop.com

Imtakt Corp./JAPAN (info@imtakt.com)

Other Countries

TI687E-KI09